labmodul



Representation of the second

LABYSTER CABINETS 60/40

Spacious with extra storage capacity


Prepared for trays and baskets

High quality for daily use

Constructed for GMO facilities

Pure design

Easy to clean

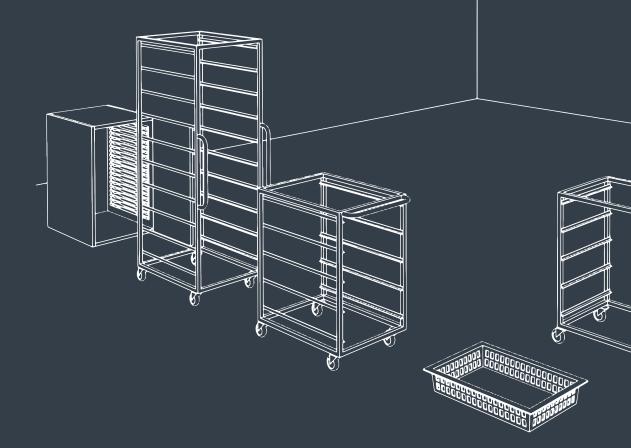



SYSTEM CABINETS 60/40

System cabinets is used for an optimum use of space inside the storage cabinets. The cabinets can be eqipped with adjustable shelves or our ISO Flex Module ® baskets and trays for maximum flexibility, efficiency.

The modular system ISO 60/40 is the result of extensive development in order to achieve maximum flexibility in systematically and safe storage that ensures order and clarity.

Universal application

System cabinets are a perfect solution for Laboratorie, hospitals, pharmacies, medical practices and wherever it is needed hygienic and mobile storage solution.

A range of mobile units for the baskets and trays are available, when safe transportation is needed.



A wide programe

System ISO 60/40 comes with a wide range of different material solutions for table tops and fronts and comes with different types of hinged doors or drawers, to fit exactly to the function of each cabinet.

Fixed dividers can easily be attached in the modules with a ClipOn - function and can be adjusted for every 10 and 20 mm. For a more flexible solution, the baskets and trays can be fitted with easy moveable dividers that can be placed in all positions.

Trays and baskets can be equipped with end-stops for fixing the basket or tray in the pull-out situation. When pulled out of the cabinet it can be positioned in a horizontal or sloped position. A basket or tray I a sloped position gives a good overview, even for those in the upper part of the cabinets. It is possible to add your own Logo or name on trays and baskets with a separate text marker. These are offered separately at request.



Trays with over 1000 options for subdivision

Baskets and trays is available in three different materials

- ABS plastic (gray)
- PC plastic (transparent)
- Wire basket (plastic coated steel)

Trays and baskets with fixed pull-out.

 Can be placed both horizontally and sloped for easy overviwe of items stored.

Trays and baskets can be fitted with an 110% pull-out bearing that withstands storage up to 50 kg.

 Ability removable rails for the absolute flexible layout of the closet

Trays and baskets with pull-out

can carry 20 kg to 50 kg

hygiene

Low weightStackable

(ClipOn effect)Ability labeling, example with company logo


Smooth surfaces for optimum

Simple installation of dividers







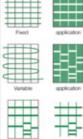

labmodul

All



Fronts

- Heighpressure laminate
- Compact laminate
- LinoliumCorian



Glass cabinet fronts

- Natural anodised 20 mm aluminum frame with 110° hinges
- White powder coated 20 mm aluminum frame with 110° hinges
- Laminate 60 mm frame with 270° hinges
- Laminate 60 mm frame concealed 110° hinges





Handles

- Standard logistics aluminum handles with inclined surface for for example Dymo strips. Width 190 mm
- Handle in brushed stainless steel.
 Width 160 mm
- Handle in matt aluminium.
 Width 160 mm







Countertops

Material options Laminate

- Compact laminate
- Stainless steel
- Corian
- Glass

Fronts



Compact laminate

Corian

Laminate colors



Laminate colors

Laminate colors

White




Yellow

Plinths

Two standard colors

Tabletop profiles

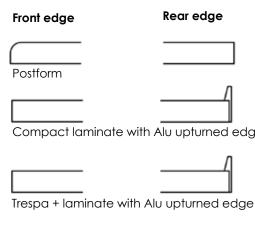



Ceramic with bead


Yellow





Compact laminate with Alu upturned edge

Terril

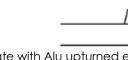


Stainless steel, with bead








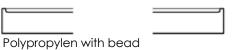


Cabinets


Low under cabinets / Widths: 456, 656 mm. Depths: 450, 650 mm. Heights 720 mm.

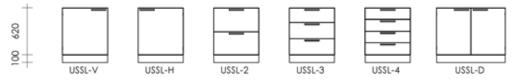

High under cabinets / Widths: 456, 650 mm. Depths: 450, 650 mm. Heights 880 mm

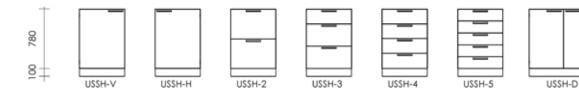


Low wall cabinets / Widths: 456, 646 mm. Depths: 350, 450 mm. Heights 600

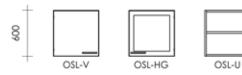






High wall cabinets / Width: 646 mm. Depth: 450 mm. Heights 750 mm

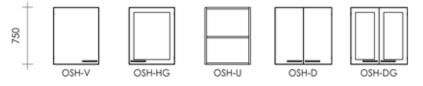


Mobile cabinets / Widths: 456, 656 mm. Depths: 450, 650 mm. Heights: 520 mm



Mobile cabinets / Widths: 456, 656 mm. Depths: 450, 650 mm. Heights: 680 mm

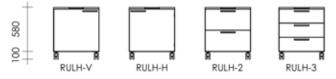


Storage cabinets / Widths: 450 mm. Depths: 450, 650 mm. Heights: 2100 mm

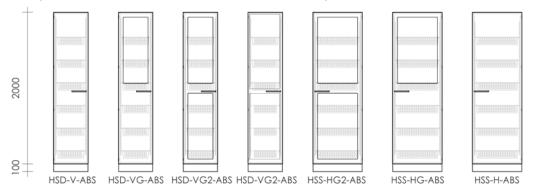



Tall storage cabinets

The individual storage system for your individual needs

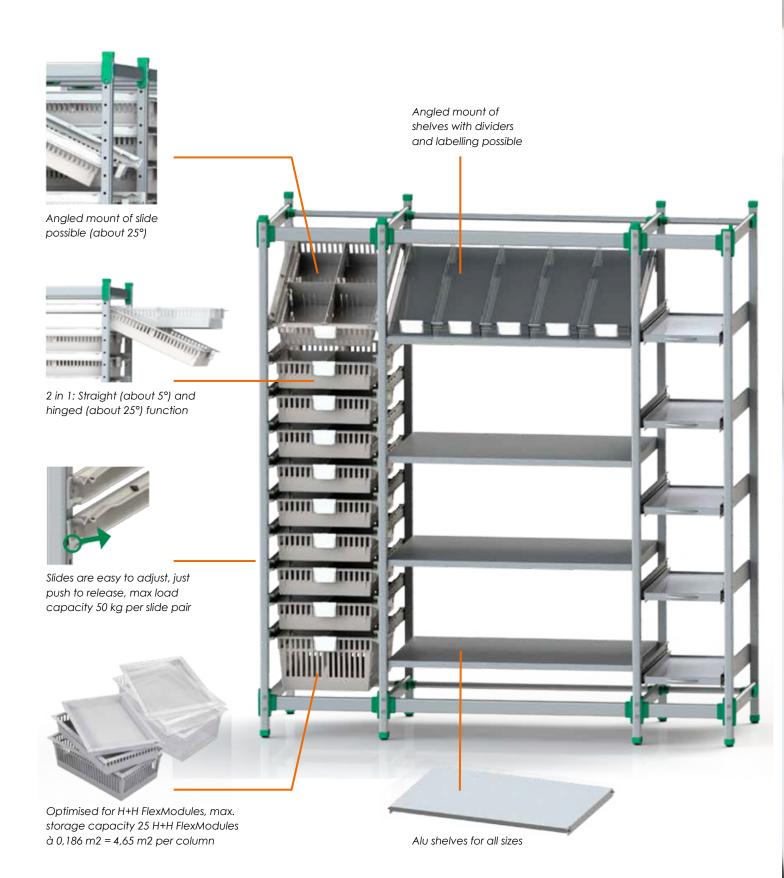



Trays and baskets

FlexModul® ISO formats stands for modular flexibility and high quality of absolute tidy storage in cabinet modules

Trays and baskets can of course be equipped with limit switches for fixed withdrawals.



See special catalog: Mobile benches and Carts Trays and Baskets

